



A total system approach to managing your breath alcohol test and calibration data.

With AlcoMark you can print, sort, edit and view your data from your computer.

The perfect tool for analyzing and archiving your results.



AlcoMark[®]


Breath Testing Management Software

Unlock the Power of Alcohol Testing

Overview

AlcoMark software was developed by Lifeloc Technologies exclusively for use with Lifeloc breath testing devices. The AlcoMark database houses breath test and calibration results from any number of Lifeloc instruments. The program enables communication between Lifeloc equipment and PCs, allowing you to view, sort, filter, print, edit and archive your breath test and calibration data. Results can be exported to Excel, XML or Text files to allow for further analysis and sharing of data.

The software is simple to install and has a user friendly interface. It is available as a single install or with a site license supporting up to five PCs.

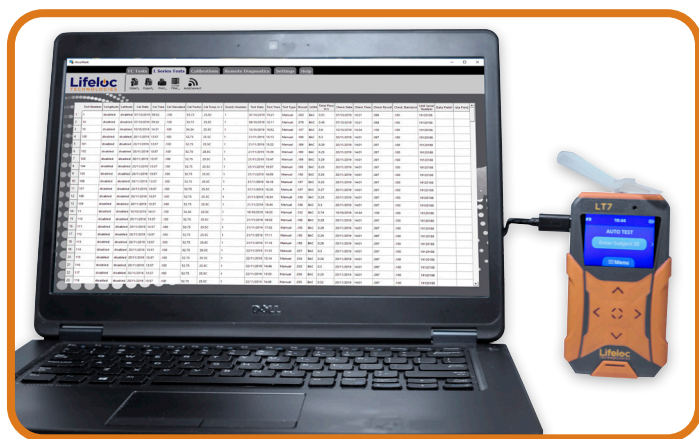

Functionality

- Load breath test and calibration results into database
- Create customized database views
- Print test and calibration results directly from AlcoMark
- Export data to Excel, XML and Text files
- Search, sort and filter results for data analysis
- Archive test and calibration results
- Import data from the EASYCAL G2 Calibration Station

Compatible Lifeloc Equipment

- LX9™
- EASYCAL® G2
- FC10^{Plus}™
- LT7™
- Phoenix 6.0BT
- FC20™
- LT71™
- EV30™
- FC20BT™

Specifications	
Software Installation	AlcoMark USB Drive
PC Operating System	Microsoft Windows 10 and 11
USB Port	USB computer cable included with AlcoMark for connecting Lifeloc device with PC
Language	English, Spanish, Portuguese
Warranty	1-Year Warranty

